

# Difficult or Emergency Airway Intubation & Ventilation Products

**New Product**

## The EasyTube



Designed for difficult tracheal or esophageal intubation. The EasyTube creates a viable airway whether it is placed in the trachea or in the esophagus. When placed in the esophagus, the EasyTube allows the use of fiber-optic devices, or passage of a suction catheter or tube exchanger. The risk of tracheal trauma is reduced to a minimum due to the single lumen at the distal tip.



- ▶ Single lumen at distal tip reduces potential for trauma
- ▶ High volume low pressure PVC cuff at distal end allows for minimal wall pressure
- ▶ Soft Latex Fee cuff to seal the pharynx
- ▶ Special placement rings clearly indicate a visible depth marking that can be used for either tracheal or esophageal intubation
- ▶ Large volume syringe, pre-filled with air for the pharyngeal cuff, and small volume syringe, pre-filled with air for the distal cuff, making intubation quick and easy


1841-41 41 Fr Large  
1828-28 28 Fr Small

The complete EasyTube set includes tube, syringes pre-filled with air, and an extra length suction catheter. Packaged in a convenient roll-up kit.

## Reusable & Disposable Endoscopy Masks

If a difficult intubation is prognoses, orotracheal or nasotracheal fiberoptic intubation becomes the chosen method. The Endoscopy Mask allows simultaneous anesthesia and ventilation of the patient. The mask can also be used for diagnostic fiberoptic Airway Endoscopy.

- ▶ The clear mask and silicone membrane allow optimal vision for introduction of the fiber bronchoscope.
- ▶ The VBM Endoscopy Mask allows mask ventilation of the patient during fiberoptic intubation.
- ▶ In connection with the VBM Endoscopy Mask, the Endotracheal Tube can be fitted with the connector, as the silicone membrane stays with the endotracheal tube after intubation, while the mask is removed over the endotracheal tube and connector.
- ▶ Diagnostic Airway Endoscopy with the fiber bronchoscope is possible using the VBM Endoscopy Mask with anesthetized patients during spontaneous breathing, assisted or controlled ventilation. The complete airway from the nose to the subsegments of the lungs can be diagnosed.
- ▶ Masks are available to aid difficult intubation for neonatal up to adult.

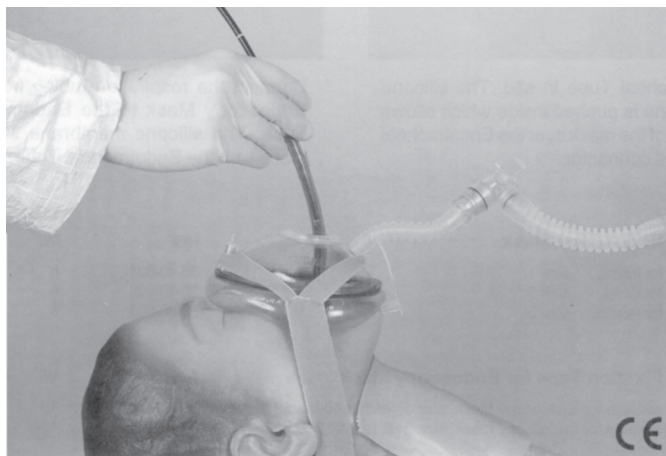

### Disposable Endoscopy Mask

6-3011	Infant	Size 1
6-3033	Child	Size 3
6-3055	Adult	Size 5

### Reusable Endoscopy Mask

6-3001	Infant	Size 1
6-3003	Child	Size 3
6-3005	Adult	Size 5

## Replacement Parts for Reusable Endoscopy Masks

### Fixation Tape for Endoscopy Mask



6-3522 Foam Tape for Mask Size 1, 400mm  
6-3533 Foam Tape for Mask Size 3 & 5


### Silicone Membrane

6-3011 Size 1 Mask  
6-3013 Size 3 Mask  
6-3015 Size 5 Mask



### Silicone Connecting Tube

15mm connector - Autoclavable up to 134°C

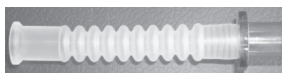

6-6098

### Mask Plug

Mask Size 1  
6-60991

Mask Size 3 & 5  
6-6099

